

Network Installation Instructions

Wavefunction Software

(Spartan Parallel Suite, Spartan Student,
Odyssey Student, Odyssey Instructor)


Installing License Server Control The license server machine can be a Windows (10 or higher, 64-bit only) or Macintosh (10.13 High Sierra or). **Use of a virtual machine as network license manager (or client machines) is available under current maintenance, annual site, or enterprise licensing. Contact sales@wavefun.com with questions**

1) Download/install/run the Wavefunction Network License Server Control Program:

[Wavefunction Downloads Page](#)


Installation will also include the current Sentinel Run-time Environment. If an older version is found you may be prompted to remove this and install the latest version.

2) Enter your Network License Product Key, and click Activate:

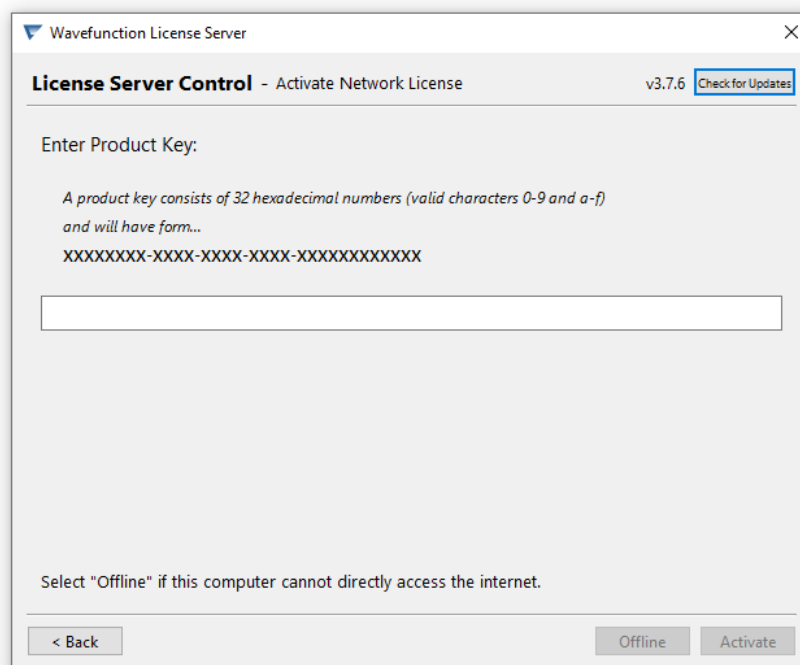

Windows version shown above

3) Successful activation will be confirmed with:

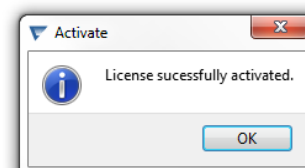

4) If you have existing Wavefunction licensing, click the Update button. You will be prompted to choose an Update method (either a Product Key or a Key Update File). Select the option you were provided with and click Next. You will either be prompted to provide your Product Key and to Activate, or to Browse and Open your Key Update file and Apply the file. Either will result in updating your licensing and the changes will be reflected in the License Server Control Program.

Configuring Client Machines: Download and install Spartan Parallel Suite, Spartan Student, or Odyssey Instructor or Student software on the client machines. Visit the [Wavefunction Downloads Page](#) and download the **MOST RECENT** versions of Spartan and Odyssey to utilize with your new licensing. For Spartan Parallel Suite customers, access [SSPD data collection downloads HERE](#).

When prompted to choose a license type, select the Network option. In the field provided, enter in the IP address of the machine that will act as the license server (see page 1). For silent installation (Windows only), see the end of this document. For silent installation under Mac OSX, email support@wavefun.com for instructions.

This completes the network installation process. Open the program to test for license access, and go to Help menu > About . . . to view license information.

Silent Installation (Windows) – For Macintosh, contact support@wavefun.com
Odyssey Instructor, Odyssey Student, Spartan Student, Spartan and Spartan Parallel Suite may be installed silently. To perform silent install:

Download and save the installer for the program you wish to install silently.
Open an Admin Level Command Prompt.
CD to the directory where the installer is saved.

For Spartan'24 Parallel Suite, run:

Spartan24Installer120.exe /S /Network=X.X.X.X
(Note: /capital letter "S" then space/)

Where X.X.X.X is the IP address of the machine that will act as the license server.

For Spartan Student, run:

SpartanStudentInstaller903.exe /S /Network=X.X.X.X

Where X.X.X.X is the IP address of the machine that will act as the license server.

For Odyssey, run:

Odyssey-7-Student-Installer.exe /S /Network=X.X.X.X —or—
Odyssey-7-Instructor-Installer.exe /S /Network=X.X.X.X

Where X.X.X.X is the IP address of the machine that will act as the license server.

Contact SUPPORT@WAVEFUN.COM with any issues.

Last Updated May 2025
wso